


Rtx[®]-TNT Columns

(6m, 0.53mm ID, 1.5µm df—cat.# 12998)

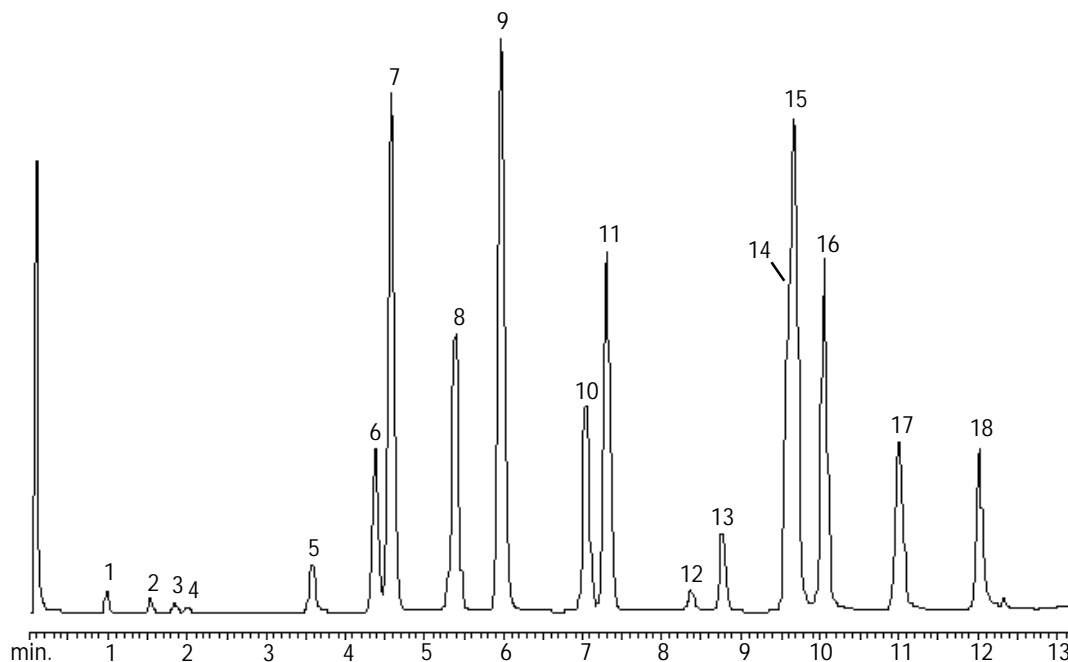

Run Conditions

On-column concentration:

200-1000pg per compound
8095 Cal. Mix A (cat.# 31607)
8095 Cal. Mix B (cat.# 31608)
3,4-dinitrotoluene (cat.# 31452)

Oven Program:

80°C (hold 1 min.) to 180°C @ 10°C/min.
to 300°C @ 30°C/min. (hold 3 min.)

Injection Port:

Direct injection using a 1mm Siltek™-
deactivated Uniliner® @ 250°C

Detector:

ECD at 330°C with anode purge

Dead Time:

4.4 seconds

Head Pressure:

He @ 3psi (20.7 KPa)

Flow Rate:

He @ 17mL/min. @ 80°C

Peak Identification

1. nitrobenzene
2. 2-nitrotoluene
3. 3-nitrotoluene
4. 4-nitrotoluene
5. nitroglycerine
6. 1,3-dinitrobenzene
7. 2,6-dinitrotoluene
8. 2,4-dinitrotoluene
9. 3,4-dinitrotoluene (IS)
10. 1,3,5-trinitrobenzene
11. trinitrotoluene
12. PETN
13. RDX
14. 4-amino-2,6-dinitrotoluene
15. 3,5-dinitroaniline
16. 2-amino-4,6-dinitrotoluene
17. tetryl
18. HMX

**Call Technical Service at 800-356-1688 or 814-353-1300, ext. 4 (or your local Restek representative)
if you have any questions about this product or any other Restek product.**



#139-01 [001]
Rev. date: 5/01

Restek USA: 110 Benner Circle, Bellefonte, PA 16823
phone: (800) 356-1688 or (814) 353-1300 • fax: (814) 353-1309; www.restekcorp.com
Restek GmbH: phone: (49) 06172 2797 0 • email: RESTEK-GMBH@t-online.de
Restek France: phone: (33) 01 60 78 32 10 • email: restekfr@club-internet.fr
Thames Restek UK Ltd.: phone: (44) 01753 624111 • email: Sales@Thamesrestek.co.uk
Thames Restek Scotland Ltd.: phone: 0870 241 1247 • email: Sales@Thamesrestek.co.uk
Restek Ireland: phone: (44)2 890 814 576 • email: dmannusrestek@compuserve.com