

Reliable Connections Made Simple


by Donna Lidgett, GC Accessories Product Marketing Manager

restek
innovation!

- Reliable seal integrity—will not unexpectedly disconnect during temperature-programmed analyses.
- Open design allows visual confirmation of the seals, for added confidence in the connections.
- Use standard Press-Tight® connectors.



Make secure, reliable column-to-column connections with SeCure™ "Y" connectors. Secondary seals ensure a leak-tight connection.

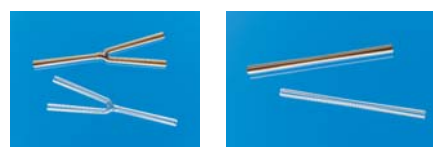

Both SeCure™ "Y" and Vu2 Union™ Connectors use standard Press-Tight® connectors—no expensive, unique inserts to purchase.

SeCure™ "Y" Connectors

Connect two analytical columns to a transfer line or guard column.

Combine the simplicity of a "Y" Press-Tight® connector with the strength of a metal union. The ferrules and knurled nuts hold the fused silica tubing in place, which prevents the tubing from unexpectedly disconnecting, even at temperatures as high as 400°C.

Kits include: SeCure™ "Y" connector body, 3 knurled nuts, 1 "Y" Universal Press-Tight® union, and 3 ferrules.

Description	Ferrules Fit Column ID	qty.	cat.#
SeCure™ "Y" Connector Kit	0.25/0.28mm	kit	20276
SeCure™ "Y" Connector Kit	0.28/0.32mm	kit	20277
SeCure™ "Y" Connector Kit	0.45/0.53mm	kit	20278
Knurled nut		3-pk.	20279

Graphite Ferrules for SeCure™ "Y" Connectors

Ferrule ID	Fits Column ID	Graphite 10-pk.	Graphite 50-pk.
0.4mm	0.25/0.28mm	20200	20227
0.5mm	0.28/0.32mm	20201	20228
0.8mm	0.45/0.53mm	20202	20224

Universal "Y" Press-Tight® Connectors

Description	ea.	3-pk.
Universal "Y" Press-Tight® Connector	20405	20406
Siltek®-treated Universal "Y" Press-Tight® Connector	20485	20486

Vu2 Union™ Connectors

Connect a guard column to an analytical column, a column to a transfer line, two columns in series, or repair a broken column.

Kits include: Vu2 Union™ body, 2 knurled nuts, 2 Press-Tight® unions, and 4 ferrules

Description	Ferrules Fit Column ID	qty.	cat.#
Vu2 Union™ Connector Kit	0.15–0.25mm	kit	21105
Vu2 Union™ Connector Kit	0.28/0.32mm	kit	21106
Vu2 Union™ Connector Kit	0.45/0.50 & 0.53mm	kit	21107
Knurled nut		2-pk.	21108

NOTE: Not recommended for GC column-to-MS connections—use the Vacuum Vu-Union® described in our catalog.



The Vu2 Union™ connector's open design allows visual confirmation of the seal; secondary seals ensure a leak-tight connection.

Graphite Ferrules for Vu2 Union™ Connectors

Ferrule ID	Fits Column ID	Graphite 2-pk.	Graphite 10-pk.
0.4mm	0.18–0.25mm	20280	20281
0.5mm	0.28/0.32mm	20282	20283
0.8mm	0.45/0.50 & 0.53mm	20284	20285

Universal Press-Tight® Connectors

Description	5-pk.	25-pk.	100-pk.
Universal Press-Tight® Connectors	20400	20401	20402
Siltek®-treated Universal Press-Tight® Connectors	20480	20449	20481