


# Extended Capillary Column Nuts

For use with standard 1/16-inch capillary ferrules - cat.# 22633, 22636

## Recommended Installation Distance:

Restek recommends an installation distance of 34 to 36 mm from the back of the nut to the end of the fused silica column.

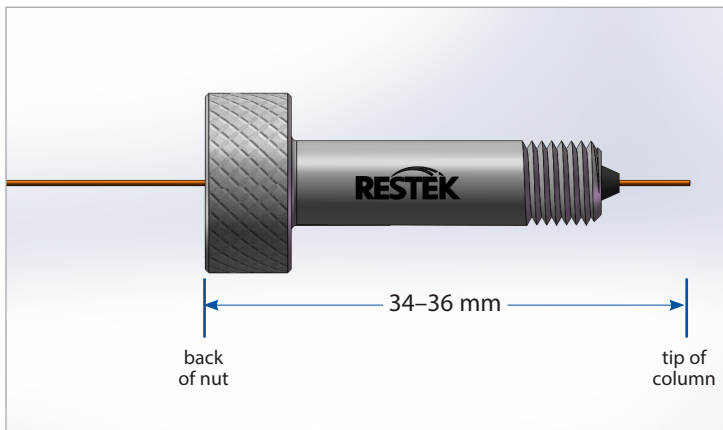

**Questions about this or any other Restek® product? Contact us or your local Restek® representative ([www.restek.com/contact-us](http://www.restek.com/contact-us)).**

Restek® patents and trademarks are the property of Restek Corporation. (See [www.restek.com/Patents-Trademarks](http://www.restek.com/Patents-Trademarks) for full list.) Other trademarks in Restek® literature or on its website are the property of their respective owners. Restek® registered trademarks are registered in the U.S. and may also be registered in other countries.

© 2014 Restek Corporation. All rights reserved.  
Printed in the U.S.A.

#203-03-040 Rev. date: 10/14



**RESTEK®**  
[www.restek.com](http://www.restek.com)