


Permanent Gases on MXT®-Msieve 5A

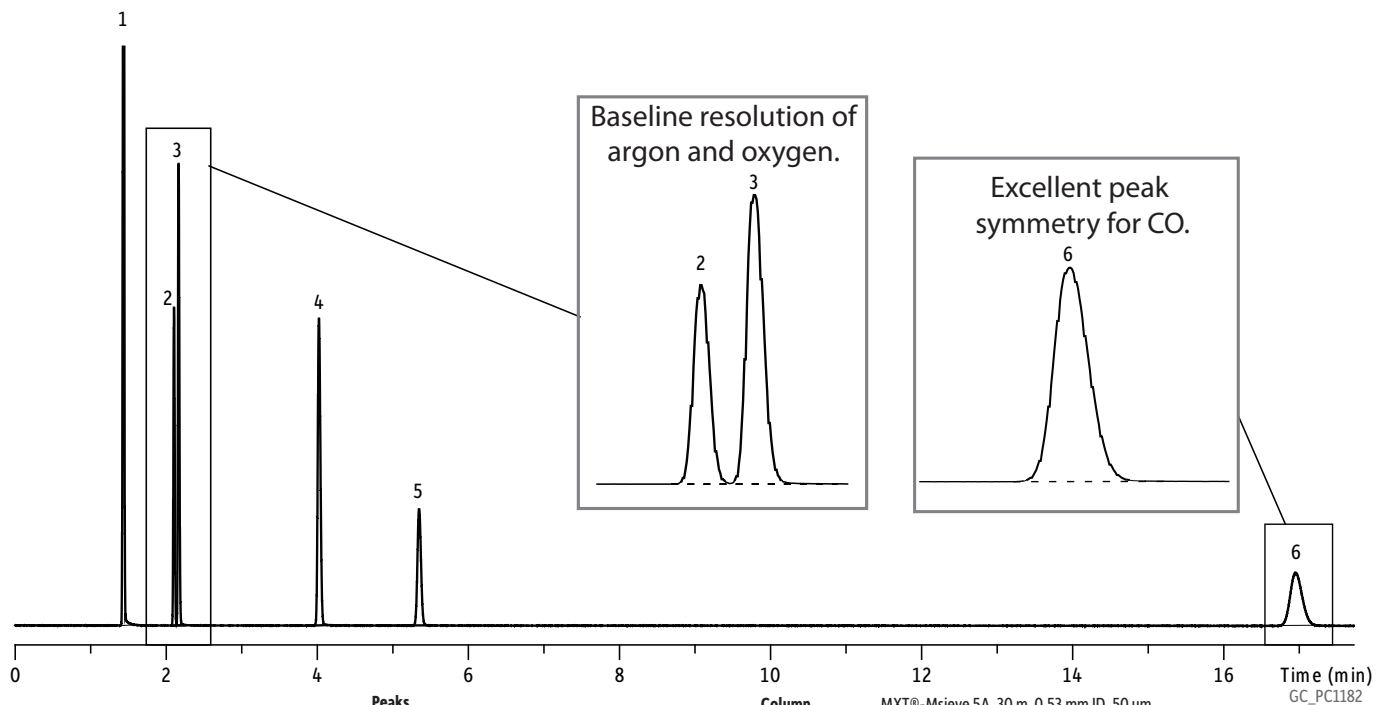

- Peaks**
1. Helium
 2. Argon
 3. Oxygen
 4. Nitrogen
 5. Methane
 6. Carbon monoxide

Column MXT®-Msieve 5A, 30 m, 0.53 mm ID, 50 µm (cat.# 79723)
Sample Conc.: 1% in hydrogen
Injection split (split ratio 50:1)
Oven
Oven Temp: 30 °C
Carrier Gas H₂
Detector µ-TCD

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
 1-814-353-1300 • 1-800-356-1688 • Fax: 1-814-353-1309 • www.restek.com