

Acrylamide Reference Standard (200 ppb) and IS on Allure Acrylamide Column

Peaks	tr (min)	Conc. (ng/mL)	Precursor	Product
1. Acrylamide-2,3,3-d3	1.739	200	75.1	58.1
2. Acrylamide	1.777	200	72.1	55.1


Column Allure Acrylamide (cat.# 9167552)
Dimensions: 50 mm x 2.1 mm ID
Particle Size: 5 µm
Pore Size: 60 Å
Guard Column: Allure Acrylamide 10 mm, 2.1 mm ID, 5 µm (cat.# 916750212)
Temp.: 22 °C
Sample Acrylamide (cat.# 30494)
Diluent: Water
Conc.: 200 ng/mL
Inj. Vol.: 10 µL

Mobile Phase

A: 0.001% Formic acid in water
B: 0.001% Formic acid in acetonitrile

Time (min)	Flow (mL/min)	%A	%B
0.00	0.4	100	0
1.00	0.4	100	0
2.00	0.4	0	100
2.01	0.4	100	0
5.50	0.4	100	0

Detector MS/MS
Ion Mode: ESI+
Mode: MRM

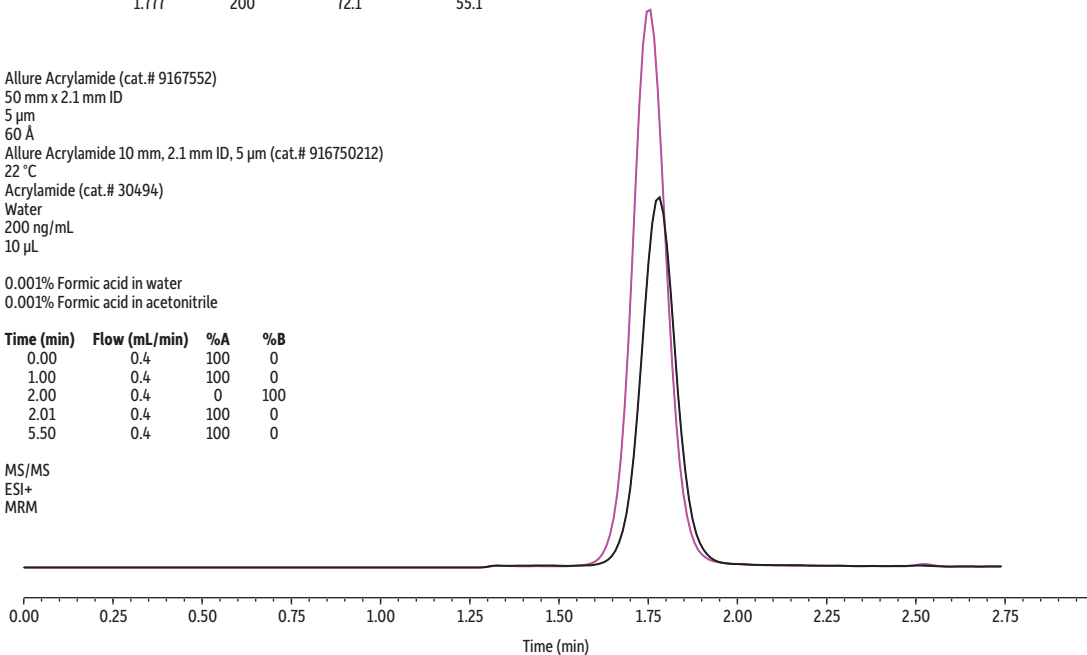

LC_FS0529



© 2019 Restek Corporation. All rights reserved.
www.restek.com

Restek patents and trademarks are the property of Restek Corporation. (See www.restek.com/Patents-Trademarks for full list.) Other trademarks in Restek literature or on its website are the property of their respective owners. Restek registered trademarks are registered in the U.S. and may also be registered in other countries.